# **Power Supply**

Model No. D2-005

Document Revision: 1

Document Last Updated on 2022/08/25 21:20

Please read Limited Warranty and General Warnings and Cautions prior to operating the D2-005.

Link to D2-005 Power Supply web page

For technical assistance, please contact Vescent at: 14998 W. 6th St. Golden, CO 80401 USA 02-005

⊟ D2-005 Power Supply

or

info [at] vescent [dot] com

# Description

The Power Supply is a quiet, linear power supply providing +/-15V and 5V. It is in a separate box to isolate 60 Hz and EM interference. One power supply can power two complete laser systems, or four electronics modules.

1/6

# Symbols

A list of symbols which appear on the D2-005, and their definitions, is provided below.



- Warning: consult manual



- Warning: risk of electrical shock. Consult manual.



- Alternating Current

– Chassis Ground

### Warnings

This manual must be consulted in all cases where the Caution symbol is marked, in order to find out the nature of the potential hazards and any actions which have to be taken to avoid them.

**Caution: Risk of Electrical Shock Do Not Open** – The D2-005 contains high voltage components and should never be opened in any circumstance except by trained Vescent personnel.



There are no serviceable parts in the D2-005. Only Vescent is authorized to service this product.



Nothing must obstruct airflow from the top panel vent. Never place anything on top of the D2-005.



A potential hazard exists if the D2-005 comes into contact with water or other conductive liquids, and exposure to these should be avoided.

Special care must be given to replacing the provided mains power cord, such that any replacement cords are appropriately rated.

If the D2-005 is used in a manner not specified by Vescent or in this manual, the protection provided by the equipment may be impaired.

# **Purchase Includes**

- D2-005 Power Supply
- VPN00464 Power Cord

# **Absolute Maximum Ratings**

Note: All modules designed to be operated in an indoor laboratory environment with filtered air. There are no known altitude restrictions on the D2-005.

Parameter	Rating
Environmental Temperature	>15°C and <30°C
Environmental Humidity	<60%

Environmental Dew Points <15°C

### **Specifications**

		Units
Current Limit		
+15V	1.5	Α
-15V	1.5	Α
+5V	6	Α
<b>Output Regulation</b>	0.05	%
Input Voltages	100±10, 120±10, 230±10	V
Input Frequency	50, 60	Hz
Power Rating	160	VA

Replaceable fuses located in the power entry module are rated at 2A, 250V @115VAC or 1A, 250V @230VAC

### Inputs, Outputs, and Controls

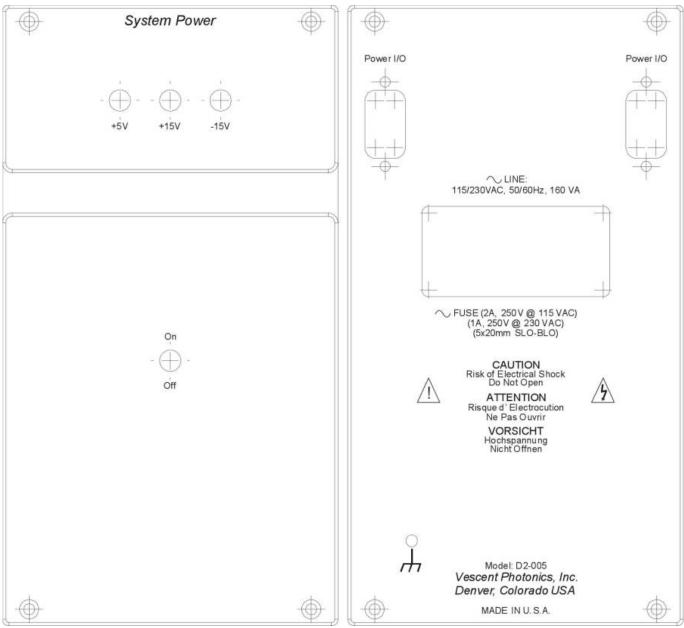


Fig. 1: Schematic drawing of the front and back panels of the D2-005

### **Monitor Section**

The monitor section contains three blue status LEDs indicating proper voltage on the +15,-15, and 5V power lines. If there is a voltage sag on any of the power lines, the blue status LED will turn off.

#### **Front Panel Controls**

A single switch on the front panel of the D2-005 toggles whether power is enabled or disabled. Flipping the switch to the "On" position turns the power supply on, enables the fan, and illuminates the three LEDs on the front. Flipping it to the "Off" position turns the power supply, fan, and LEDs off.




Do not position the D2-005 such that it is difficult to operate the switch.

#### Back panel I/O

#### Power I/O (9-pin D-sub)

Always turn off the D2-005, make connections to the device, and then re-energize the D2-005. Never connect this D2-005 Power Supply to a device when it is switched on and supplying power.

Power is supplied to the D2-005 via a standard C13 to Type-B power cord. The power to each electronics module is through a 9-pin D-sub connector through any serial cable with 9-pin D-sub connectors, which is convenient when the unit must be taken out of line for access to the side panels. There are two 9-pin D-sub connectors on the back of the D2-005. Either or both can be used to power external electronics, however, current is shared between the two connectors. The maximum power rating for the 5V line is 30W (6A), and for the  $\pm$ 15V line is 22.5W (1.5A). The pin outs are shown in the following figure:

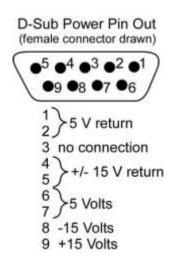

In some system configurations, the D2-5 can radiate observable line noise on nearby D2 instruments and other electronics. This will appear as a signal at the frequency of your mains electricity (typically either 50 Hz or 60 Hz). This noise can easily be removed by moving the D2-005 at least 18 inches (45cm) away from other electronics, rotating it 90° such that the sides of the D2-005 face away, or by moving the entire power supply to a different shelf. To accommodate this, all D2-005's are shipped with a 5' DB9 cable as of January 1, 2022.

If you are considering using an Uninterruptible Power Supply (UPS) to guarantee operation of your Vescent products through a brown out or black out, great care should be taken in choosing the model. Lower cost models tend to produce modified square wave voltage profiles. The high-frequency components of such a voltage profile may interact poorly with the D2-005 (and down-stream active Vescent modules). If you do choose to use a UPS, select a model that will



provide a sine wave voltage profile without higher harmonic components to avoid potential damage to your high-value equipment.

From: https://www.vescent.com/manuals/ - **Product Manuals** 

Permanent link: https://www.vescent.com/manuals/doku.php?id=d2:power\_supply&rev=1661462446



Last update: 2022/08/25 21:20